

■ ■ wissner-  
■ ■ bosserhoff

Functional elegance in nursing care

SafeSense®

Bed exit assistance system



# Bed exit assistance system

What exactly does this mean?


SafeSense® signals the nursing staff that the resident/patient has left the bed or is about to leave it. This way, falls can be detected early on and their consequences can be reduced. Essentially, the bed-exit assistance system works like a traffic light. It is very easy to understand:



Timer not activated



Timer activated  
→ e.g. 15 minutes to go to the restroom at night

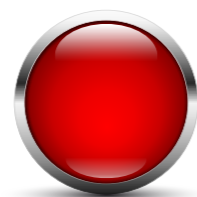

Timer elapsed  
→ Signal to the nurse call system

# SafeSense® = Safety + freedom

By day and by night – thanks to the intelligent bed exit sensor system

Sustainable fall prevention for utmost safety. This is what the SafeSense® bed exit assistance system of wissner-bosserhoff stands for.

On a national average, the annual number of falls in nursing homes and hospitals is 1.5 falls per bed\*. The consequences of falls to a patient's health are often very serious for elderly people and, in the worst case, even irreversible.

SafeSense® can considerably increase the level of safety, especially during night and weekend shifts with fewer staff. Especially at night, the bed exit assistant system can minimize the need for side guards. The freedom of the resident/patient is upheld. An automatic night light helps the resident/patient to safely find the way to the bathroom at night. It is automatically switched on when the resident/patient leaves the bed and is switched off on return. The interval of the nurse call timer can be set flexibly according to the individual degrees of mobility of the residents/patients.

**5**  
Million cases of elderly people falling each year:

Almost one third of people aged 65 and above, as well as 50% of those aged 80 and above, fall at least once a year.  
In Germany alone, there are about 5 million cases of falls among the elderly population each year.

(Source: "Gesundheit und Krankheit im Alter" (health and illness in old age). Publisher: Karin Böhm, Federal Statistical Office of Germany / Clemens Tesch-Römer, German Center of Gerontology / Thomas Ziese, Robert Koch Institute, Berlin 2009)

↓

Prevention with sentida low nursing beds

↓

Prevention with image 3 for universal hospital bed with low position

# Fall prevention in nursing homes and hospitals



## Greater safety – fewer falls/restraints

- Faster response time in the event of falls
- Greater safety during night and weekend shifts
- Ideal in combination with sentida and image 3 low nursing beds and SafeFree® side guards
- Systematic implementation of the findings from the initiative for the reduction of restraints, the "Werdenfeler Weg" (Werdenfels care pathway)



## Preserving mobility

- By using low nursing beds "sentida" / "image 3"
- With the WIBO 3-stop strategy
- With WIBO Safe-Free® side guards



## Individual protection

- At night, the way to the bathroom is lit with a sensor-controlled night light
- The nurse call timing is flexibly set from 0 to 30 min. according to the mobility of the resident/patient



## The advantages for your institution

- Investment security
- Patient/resident safety in any working position
- Legal certainty without having to comply with additional requirements
- Works effectively and efficiently
- No additional costs for safety equipment
- Efficient use of restraining belts



# How to order SafeSense®

For WIBO nursing home beds



- Integrated as an option
- As retrofit kit for all electrical nursing home beds by wissner-bosserhoff (manufactured in 2002 or later)
- Available with cable or as a wireless version


For hospital bed image 3



- Can be ordered as an option
- Available/recommended with cable or as a wireless version

## SafeSense® components


Sensor mat 2

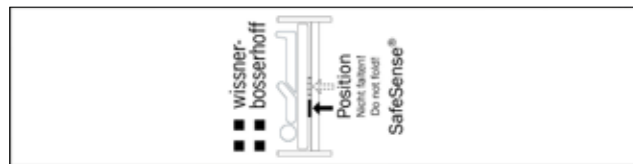


02-000781  
80 x 20 cm

Connection kit


Adapter cable with and without on/off switch and compatible plug to connect it to the nurse call system (customer-specific), for cable and wireless version

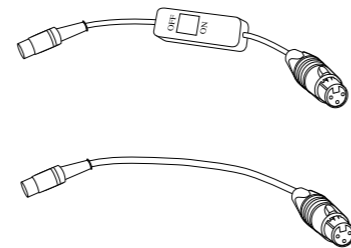

Bed-Exit-Box 2

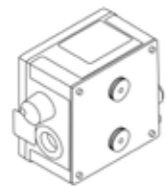

04-1851  
with universal holder  
Available with cable or as a wireless version

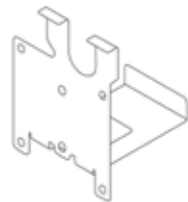

01-011193  
Retaining plate  
for DU09 motor



01-011194  
Retaining plate  
for Q7 motor

# Wireless receiver

For nursing home beds only!

 <p>Call receiver for nurse call systems</p> <ul style="list-style-type: none"> <li>• Compatible with all standard nurse call systems using the WIBO plug adapter</li> <li>• Mini receiver</li> <li>• 32 transmitters can be registered → individual room furnishing according to the resident's wishes/requirements</li> <li>• With nursing care function. Call is locked for 15 min.</li> </ul>	 <p>Plug pager</p> <ul style="list-style-type: none"> <li>• Acoustic and visual nurse call signals</li> <li>• Easy acknowledgment directly on the plug pager</li> <li>• Perfect for use in smaller residential homes and groups without a nurse call system</li> <li>• Energy-efficient! Needs only 0.3 W in stand-by mode</li> </ul>	 <p>Plug receiver</p> <ul style="list-style-type: none"> <li>• Switches on the orientation light or floor lamp in the room, for example</li> <li>• 32 transmitters can be registered</li> <li>• Sensor signal is automatically transmitted on bed exit</li> <li>• Energy-efficient! Needs only 0.2 W in stand-by</li> </ul>	 <p>Alarm telephone including wristband transmitter</p> <ul style="list-style-type: none"> <li>• Home emergency call center and comfortable telephone in one</li> <li>• Received call is successively forwarded to three pre-programmed numbers</li> <li>• Talk to the resident for three minutes with the integrated hands-free function</li> <li>• Illuminated display</li> </ul>
--	--	--	--

## Functional principle



## Application recommendations

		RH* / NH* <u>with</u> nurse call system	Small RHs and NHs and residential groups <u>without</u> nurse call system	Home care
Bed exit box 2	Cable	✓		
	Wireless	✓	✓	✓
Call receiver		✓		
Plug pager			✓	✓
Plug receiver		✓	✓	✓
Alarm telephone including wristband transmitter			✓	✓

Members of LINET Group

Get in contact:



wissner-bosserhoff GmbH

Hauptstraße 4-6 | 58739 Wickede (Ruhr) | Germany

Phone: +49(0)2377.784-0 | Fax: -167 | info@wi-bo.de | www.wi-bo.com



EU patent no.:  
EP000001623654  
EP000002052707

